


phi2 power balance — time-domain normalized scope traces — 215818_power-balance_final-bs-inputs-long

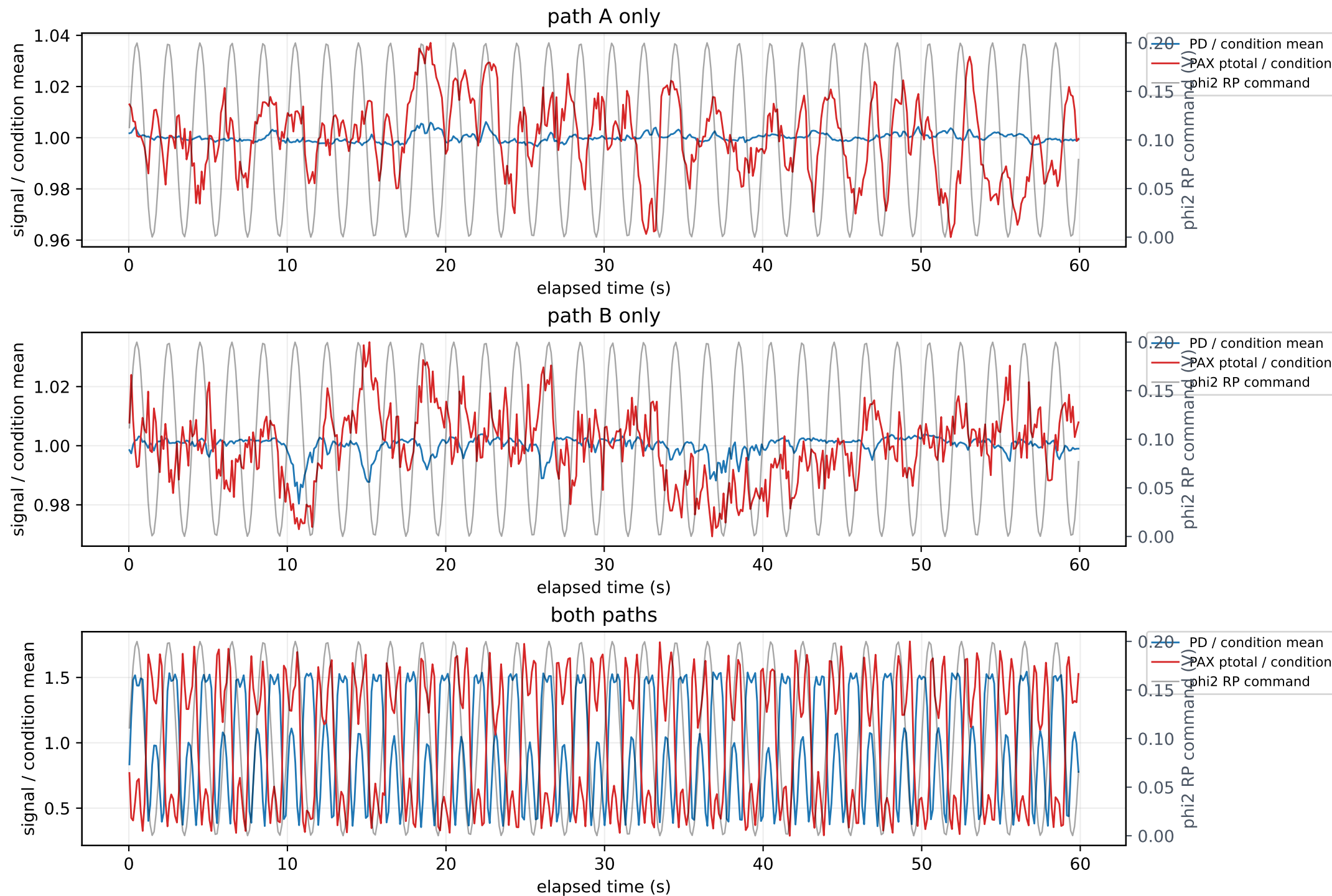

Phi2 power-balance calculations — 215818_power-balance_final-bs-inputs-long

condition	samples	mean PD (V)	PD σ (V)	PD contrast (5-95%)	mean PAX ptotal	PAX σ	PAX contrast (5-95%)
path A only	499	0.0958621	0.000156	0.0027	0.000126761	1.9e-06	0.0246
path B only	498	0.110823	0.000353	0.0047	0.00010684	1.29e-06	0.0208
both paths	498	0.230093	0.0934	0.6030	0.000229513	0.00011	0.6500

Contrast uses robust 5th/95th-percentile extrema. PD volts and PAX ptotal are different detector units; a physical 40/60 port split requires a relative detector-power calibration.